


# ETHOS™ LP

## LOW PROFILE IN HIGH GEAR



MOVEMENT  
MADE  
EASY

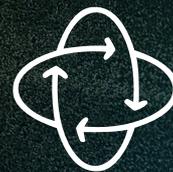

#HUMANFIRST



# ETHOS™ LP

## FEATHERWEIGHT CHAMPION MASTER OF COMPLIANCE

**Low profile. Low weight. High compliance.**  
ETHOS LP combines an exceptionally **low build height** with excellent terrain **adaptability**, natural **flexibility** and great multiaxial compliance. Delivering a **smooth, reliable and confident** gait on any surface.

Designed for users who prioritize natural movement, all-day comfort and wearability on any surface.



**Half rocker sole plate:** Designed to provide constant contact with the ground for a smooth rollover and superior ground compliance



**Split keel and heel:** Adapts to uneven terrain, slopes, and objects on the ground



**Waterproof:** Resistant to fresh, sea, and chlorinated water



**No bolts** eliminates dead spots.



**Heel bumper** provides shock absorption at heel strike.



### FUSIONFLEX COMPOSITE TECHNOLOGY

A lightweight three-layer (Oreo) composite design that balances:  
- Durability and flexibility with a fiberglass core  
- Responsiveness with carbon fiber outer layers



**MOVEMENT MADE EASY**  
± 8° of axial rotation  
± 16.5° of inversion and eversion

### MULTIAXIAL ROTATION PLATE

A proprietary rubber blend specifically designed for increased comfort, confidence, and freedom of movement used in an independent feature that allows axial rotation, inversion, eversion, and separation of the heel and keel. Designed to provide excellent flexibility for **greater ground compliance**, and **reduction of internal and external rotational forces**.



# ETHOS™ LP

LOW PROFILE IN HIGH GEAR

## ADULTS



Available in  
**EVAQ8®**

Item N°.	Size(cm)	Build height	Weight*	Categories	Side	Sandal toe	Heel height
LP3-00	22 to 25	57 mm	425 g	1 to 9	R / L	yes	10 mm
	26 to 28	58 mm	485 g				
	29 to 30	63 mm	656 g				
LP3-V3	26 to 28	58 mm	525 g				
	29 to 30	63 mm	696 g				



\*26cm Cat 4, module with foot shell

Categories based on patient's weight and impact level for ETHOS LP for adults

Weight/ Impact level	44-52kg	53-59kg	60-68kg	69-77kg	78-88kg	89-100kg	101-116kg	117-130kg	131-147kg	148-166kg
Low	1	1	2	3	4	5	6	7	8	9
Moderate	1	2	3	4	5	6	7	8	9	-
High	2	3	4	5	6	7	8	9	-	-

LP3-\_\_O\_\_A\_\_S\_\_  
**Model**  
 Standard : 00  
 EVAQ8 : V3  
**Category** 1-9  
**Size** 22-30  
**Side**  
 R- right  
 L- left

FTC-2K-1\_\_4-S\_\_  
**Size** 22-30  
**Skin tone**  
 L-Light  
 M-Medium  
 D-Dark  
**Side**  
 R- right  
 L- left



## PEDIATRICS



Item N°.	Categories	Size (cm)	Sandal Toe	Heel height	Build height	Weight
LP3-00	1P to 9P	19 to 21	no	10 mm	42 mm	275 g

Categories based on patient's weight for ETHOS LP for pediatrics

Weight/ Impact level	<26kg	26-36kg	36-45kg	46-60kg
High	1P	2P	3P	4P

\*20 cat. 2P, with foot shell

LP3-00-0\_\_P\_\_RU  
**Category** 1-4  
**Size** 19-21

FTC-2M-1\_\_4-R\_\_  
**Size** 19-21  
**Skin tone**  
 L-Light  
 M-Medium  
 D-Dark  
**Side**  
 R- right  
 L- left

